



FitMax™ CPAP Total Face Mask

for use with CPAP and bilevel NIV treatment




Critical Care ▪ FitMax™ CPAP Total Face Mask

FitMax™ CPAP Total Face Mask

The FitMax CPAP total face mask is an adult (>30kg) single patient use NIV patient interface for use with CPAP and bilevel treatment within the hospital environment. The design of the mask helps to eliminate the nasal bridge issues that can be associated with traditional oronasal masks. The FitMax is easy to fit and is available in 3 sizes for the majority of patients.

Two colour coded elbow options are provided for different clinical situations – an anti-asphyxia elbow with leak port for use with single limb breathing systems and a standard elbow for use with dual limb breathing systems.

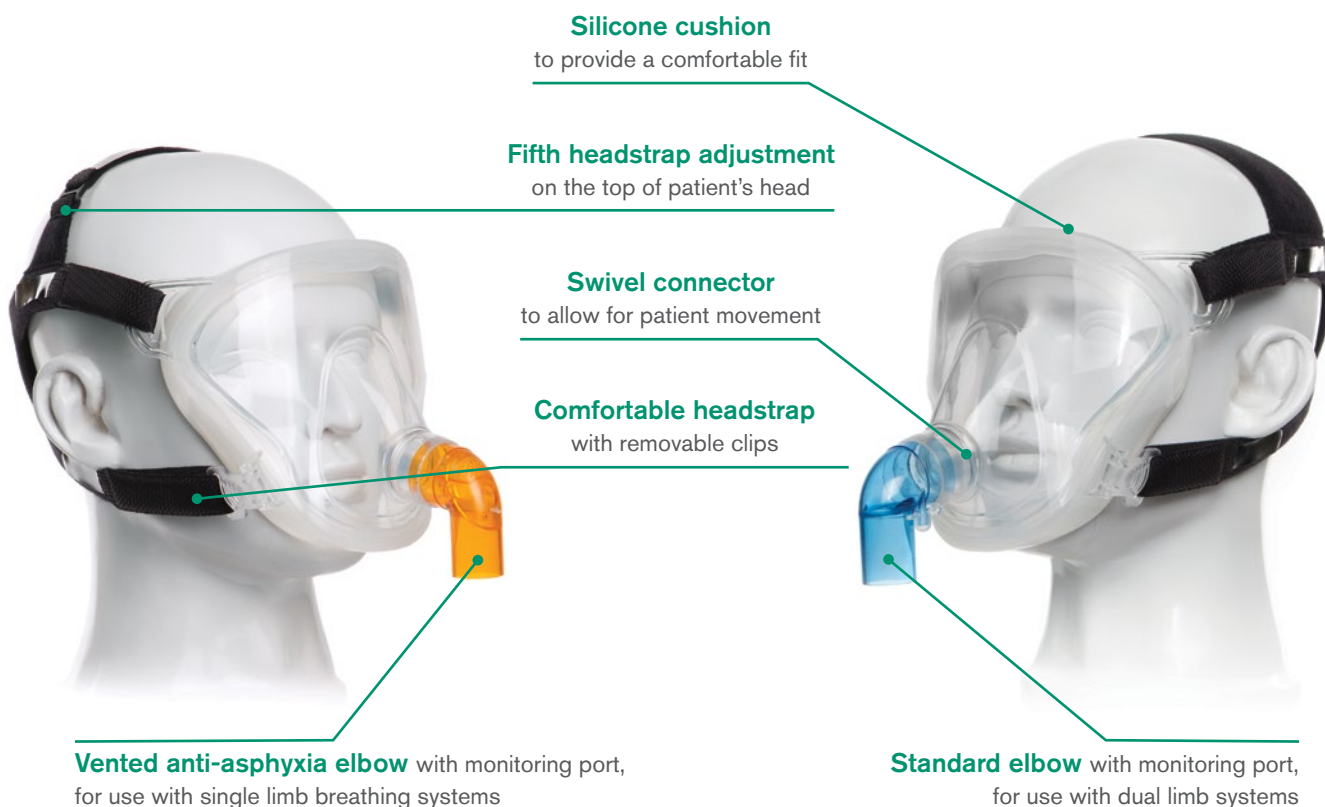

Code	Description	Colour coded elbow	Box Qty.
For use with single limb breathing systems			
2290020	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, medium	●	10
2290030	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, large	●	10
2290040	FitMax CPAP total face mask, with vented anti-asphyxia elbow, head strap, extra large	●	5
For use with dual limb breathing systems			
2290070	FitMax CPAP total face mask, with standard elbow, head strap, medium	●	10
2290080	FitMax CPAP total face mask, with standard elbow, head strap, large	●	10
2290090	FitMax CPAP total face mask, with standard elbow, head strap, extra large	●	5

[Make an enquiry](#)

The FitMax masks are manufactured by Besmed Health Business Corporation, Taiwan and distributed by Intersurgical Ltd

Please think before you print Save energy and paper. If you must print this information sheet please print it double sided.

IS1.52 • Issue 1 01.21



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com

UK • Ireland • France • Germany • Spain • Portugal • Italy • Benelux • Sweden • Denmark • Lithuania • Russia • Czech Republic
Turkey • South Africa • China • Japan • Taiwan • Philippines • USA • Canada • Colombia • Australia