



GO-PAP™

Disposable CPAP device



Critical Care ▪ Pulmodyne®


GO-PAP™

The GO-PAP is a disposable emergency CPAP device that connects directly to an oxygen source. It operates at 10 L/min and delivers approximately 30% FiO₂. There is an adjustable PEEP valve (5, 7.5 and 10 cmH₂O) integrated into the generator. The PEEP is independent of the flow and therefore

there is no requirement to adjust the flow to maintain PEEP levels. Additionally, this completely disposable system offers integrated nebulisation while maintaining a consistent flow to the patient.

Two different headgear designs

standard or CAP headgear (as shown)

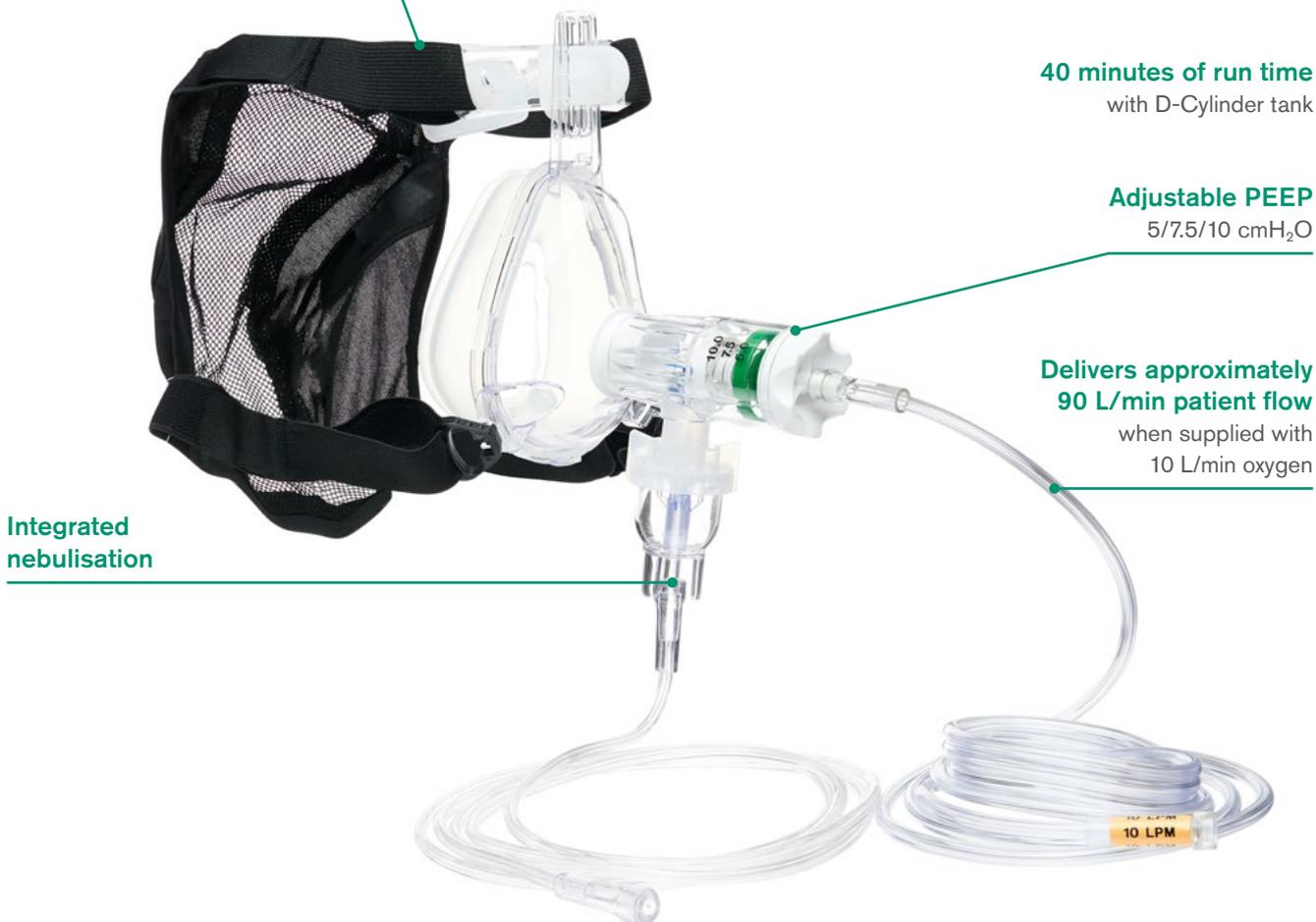

40 minutes of run time
with D-Cylinder tank

Adjustable PEEP
5/7.5/10 cmH₂O

Delivers approximately 90 L/min patient flow
when supplied with 10 L/min oxygen

Integrated nebulisation

GO-PAP System with nebuliser and BiTrac™ ED mask with CAP headgear

Code	Description	Box Qty.
313-4605N	GO-PAP System with oxygen tubing, nebuliser and BiTrac ED mask with CAP headgear, large adult	10
313-4604N	GO-PAP System with oxygen tubing, nebuliser and BiTrac ED mask with CAP headgear, medium adult	10



GO-PAP™ System and BiTrac™ ED mask with CAP headgear

Code	Description	Box Qty.
313-4605	GO-PAP System with oxygen tubing and BiTrac ED mask with CAP headgear, large adult	10
313-4604	GO-PAP System with oxygen tubing and BiTrac ED mask with CAP headgear, medium adult	10



GO-PAP System with nebuliser and BiTrac ED mask with standard headgear

Code	Description	Box Qty.
313-4603N	GO-PAP System with oxygen tubing, nebuliser and BiTrac ED mask, large adult	10
313-4602N	GO-PAP System with oxygen tubing, nebuliser and BiTrac ED mask, medium adult	10
313-4601N	GO-PAP System with oxygen tubing, nebuliser and BiTrac ED mask, small adult	10



GO-PAP System and BiTrac ED mask with standard headgear

Code	Description	Box Qty.
313-4603	GO-PAP System with oxygen tubing and BiTrac ED mask, large adult	10
313-4602	GO-PAP System with oxygen tubing and BiTrac ED mask, medium adult	10
313-4601	GO-PAP System with oxygen tubing and BiTrac ED mask, small adult	10



BiTrac™ ED mask with CAP headgear

Code	Description	Box Qty.
313-7031C	BiTrac ED mask with CAP headgear, extra large adult	10
313-7030C	BiTrac ED mask with CAP headgear, large adult	10
313-7029C	BiTrac ED mask with CAP headgear, medium adult	10
313-7028C	BiTrac ED mask with CAP headgear, small adult	10



BiTrac ED mask with standard headgear

Code	Description	Box Qty.
313-7031	BiTrac ED mask, extra large adult	10
313-7030	BiTrac ED mask, large adult	10
313-7029	BiTrac ED mask, medium adult	10
313-7028	BiTrac ED mask, small adult	10

GO-PAP™ accessories

Accessories include the GO-PAP System generator and the ADAPT-2-VENT™ which is an adapter elbow that allows the use of the BiTrac ED mask on standard dual limb or single limb ventilator circuits. This allows the patient to move from CPAP, using the GO-PAP™ or the O2-MAX™, to other forms of ventilation without changing the patient interface.



313-4600



313-9185



313-9186

Code	Description	Box Qty.
313-4600	GO-PAP System generator with oxygen tubing	10
313-9185	ADAPT-2-VENT 30mmF non-vented elbow with anti-asphyxiation valve	10
313-9186	ADAPT-2-VENT 30mmF standard non-vented elbow	10

ADAPT-2-VENT, BiTrac, GO-PAP and O2-MAX are manufactured by Pulmodyne Inc, USA and distributed by Intersurgical Ltd

Please think before you print Save energy and paper. If you must print this information sheet please print it double sided.

IS1.63 ▪ Issue 2 07.22



Intersurgical Australia Pty Ltd, 4/151 Beauchamp Road, Matraville, NSW, 2036, Australia
 T: +61 (0)2 8048 3300 info@intersurgical.com.au www.intersurgical.com.au