



StarMed CaStar R Up

Openable hood for non-invasive mechanical ventilation (NIV)



Critical Care ▪ Patient Interfaces



StarMed CaStar R Up openable hood for NIV

Openable CaStar R Up is a comfortable, versatile and lightweight patient interface for non-invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R Up is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for NIV therapy.


Advantages

- Allows rapid access to the patient during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood

Technical Features

- 22M connectors compliant to ISO standard
- Built in inflatable neck cushion for added patient comfort
- Comfortable and atraumatic underarm fastening straps
- Other fastening accessories are available if required
- Two new larger access ports on the collar for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Read more evidence on www.intersurgical.com.au

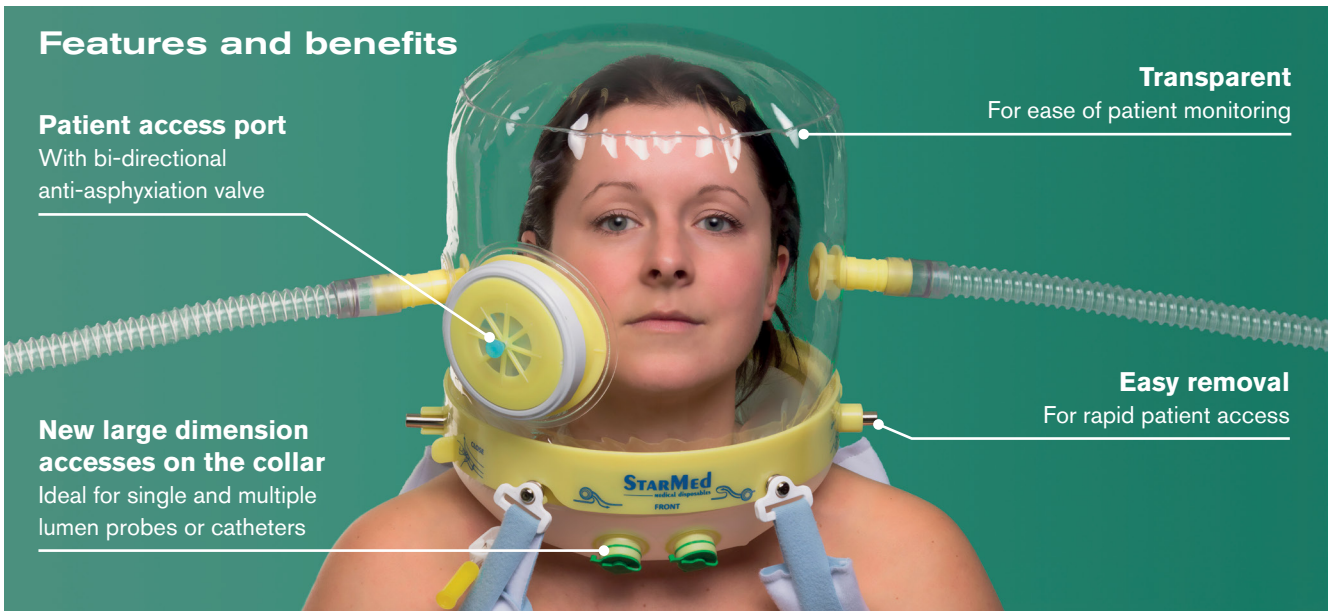

CaStar R Up openable hood for NIV with bi-directional anti-asphyxiation valve and patient access port, inflatable neck cushion and inflation bulb, ear plugs and measuring tape

Code	Size	Neck circumference	Box Qty.
CA022XXL/2RW	Extra Extra Large	≥50cm	5
CA022XL/2RW	Extra Large	45-52cm	5
CA022L/2RW	Large	40-47cm	5
CA022M/2RW	Medium	34-41cm	5
CA022S/2RW	Small	27-34cm	5
CA022XS/2RW	Extra Small	17-27cm	5



www.intersurgical.com.au/info/starmed

Make an enquiry

IS1.39_AU • Issue 3 06.16



Unit 4, 151 Beauchamp Road, Matraville, NSW, 2036

T: +61 (0) 2 8048 3300 info@intersurgical.com.au www.intersurgical.com.au

UK • Ireland • France • Germany • Spain • Portugal • Italy
 Netherlands • Sweden • Denmark • Lithuania • Russia
 Czech Republic • Turkey • South Africa • China • Japan
 Taiwan • Philippines • USA • Colombia • Australia



Intersurgical Ltd is certified to ISO 14001, ISO 9001 and ISO 13485

Please think before you print
 Save energy and paper.
 If you must print this e-mail
 please print it double sided