


# StarMed CaStar Up

Openable hood for CPAP therapy



Critical Care ▪ Patient Interfaces

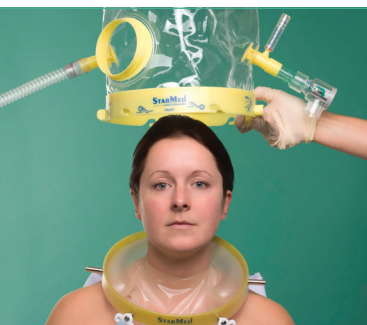

## StarMed CaStar Up openable hood for CPAP therapy

Openable CaStar Up is a comfortable, versatile and lightweight CPAP patient interface designed to offer an alternative to traditional methods of CPAP delivery. The CaStar Up is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for CPAP therapy.

### Advantages

- Allows rapid access to the patient during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood

### Technical Features

- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Built in inflatable neck cushion for added patient comfort
- Comfortable and atraumatic underarm fastening straps
- Other fastening accessories are available if required
- Two new larger access ports on the collar for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Read more evidence on [www.intersurgical.com.au](http://www.intersurgical.com.au)

**Features and benefits**

**Patient access port**  
With bi-directional anti-asphyxiation valve

**New large dimension accesses on the collar**  
Ideal for single and multiple lumen probes or catheters

**Transparent**  
For ease of patient monitoring

**Integrated pressure manometer**

**Easy removal**  
For rapid patient access

STARMed  
COMPLETE RESPIRATORY SYSTEMS

CaStar Up openable CPAP hood with bi-directional anti-asphyxiation valve and patient access port, inflatable neck cushion and inflation bulb, ear plugs and measuring tape

Code	Size	Neck circumference	Box Qty.
CA012XXL/2W	Extra Extra Large	≥50cm	5
CA012XL/2W	Extra Large	45-52cm	5
CA012L/2W	Large	40-47cm	5
CA012M/2W	Medium	34-41cm	5
CA012S/2W	Small	27-34cm	5
CA012XS/2W	Extra Small	17-27cm	5



[www.intersurgical.com.au/info/starmed](http://www.intersurgical.com.au/info/starmed) Video available at [www.intersurgical.com.au](http://www.intersurgical.com.au)

Make an enquiry

IS1.31\_AU • Issue 4 06.16



Unit 4, 151 Beauchamp Road, Matraville, NSW, 2036

T: +61 (0) 2 8048 3300 [info@intersurgical.com.au](mailto:info@intersurgical.com.au) [www.intersurgical.com.au](http://www.intersurgical.com.au)

UK ▪ Ireland ▪ France ▪ Germany ▪ Spain ▪ Portugal ▪ Italy  
Netherlands ▪ Sweden ▪ Denmark ▪ Lithuania ▪ Russia  
Czech Republic ▪ Turkey ▪ South Africa ▪ China ▪ Japan  
Taiwan ▪ Philippines ▪ USA ▪ Colombia ▪ Australia



**Please think before you print**  
Save energy and paper.  
If you must print this e-mail  
please print it double sided

Intersurgical Ltd is certified to ISO 14001, ISO 9001 and ISO 13485